


The QUEST Table Top Cassette Holder (Q-CHT) is a great accessory to any X-Ray system. It is designed to accommodate portable DR Panels, CR, and Film Cassettes ranges from 7" – 17" (18cm – 43cm), with convenient size adjustment locks which allows for the imaging plates to be positioned either portrait or landscape.

The Table Top Cassette Holder is simply positioned on the table top surface. The base is lined with felt material to protect the table top surface from scratches. The light weight design also provides the ability to optimize positioning quickly and effortlessly while providing an excellent platform for assisting in acquisition of wide range of radiographs.


From Emergency Rooms to Radiographic rooms in Hospitals, Orthopedic Facilities, Clinics and the Private practice, the Table Top Cassette Holder can assist and reduce positioning time which can be critical in certain situations, giving overall better examination results.



Shown with Wireless DR Panel
In Landscape Position

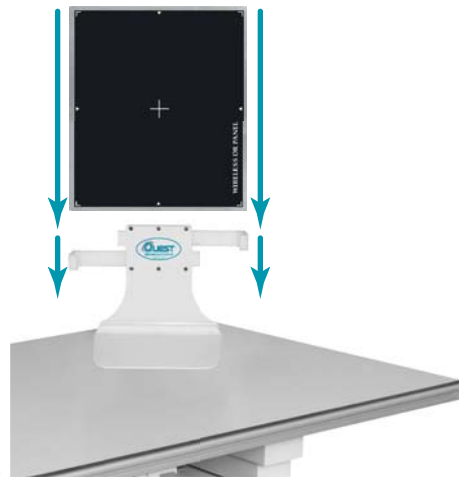

Shown with Wireless DR Panel
In Portrait Position

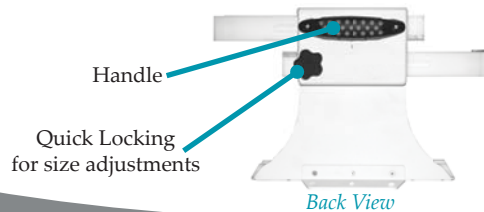

Quick Locking
for size adjustments

Back View

KEY FEATURES

CAT# Q-CHT

TABLE TOP CASSETTE HOLDER

- Quick and Easy positioning on the X-Ray table
- Accommodates all portable digital (DR) panels up to 17" x 17" (43cm x 43cm)
- Accommodates all CR and/or Film Cassettes up to 17" x 17" (43cm x 43cm)
- Quick locking side adjustments hug the imaging plate in place allowing for the imaging plates to be positioned either horizontally or vertically
- Assists in acquiring cross-table lateral exams, as well as many other specialty views
- Effortless maneuverability with table top scratch protection

*Specifications are subject to change



Improve Your Positioning For
Faster & Easier Examinations

www.questxray.com



CASSETTE HOLDERS

For Portable DR Panels as well as all CR and Film Cassettes

The QUEST Table Mounted Cassette Holder (Q-CHTM) is specifically designed to be used with portable Digital panels, as well as all CR and film cassettes up to 17" x 17" and easily slides onto the side rails of many radiographic tables. *Consult Factory

The Table Mounted Cassette Holder's articulating swivel arm allows the technologist to position it close to the patient for cross-table, as well as decubitus exams. The $\pm 90^\circ$ rotation of the articulating swivel arm allows for angle views of various anatomical regions. The technologist no longer needs to stand in the exam room to hold the plate in position.




KEY FEATURES

CAT# Q-CHTM

TABLE MOUNTED CASSETTE HOLDER

- Articulating swivel arm ($\pm 90^\circ$) allows variety of angled positions, along table-top
- Accommodates portable DR panels, CR or film cassettes, up to 17" X 17"
- Mounts to either side of Table Rails

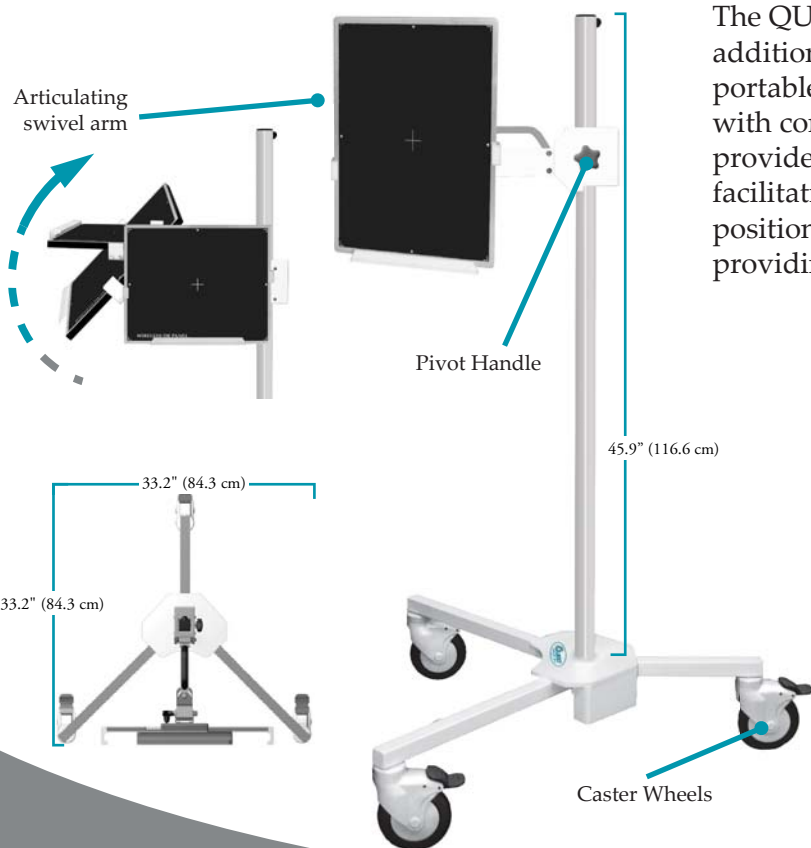

The QUEST Rolling Cassette Holder (Q-CHR) is another great addition to any X-Ray examination facility. It accommodates all portable DR Panels, CR and Film Cassettes up to 17" x 17" in size with convenient size adjustment locks. The ergonomic design provides the ability to position imaging plates around the patient facilitating a precise and stable imaging platform & reduces positioning time which can be critical in some areas, along with providing an excellent platform for performing patellar views.

KEY FEATURES

CAT# Q-CHR

ROLLING CASSETTE HOLDER

- Accommodates portable DR panels, CR or film cassettes, up to 17" X 17"
- Quick and Secure maneuverability around the X-Ray exam room
- Height Adjustable positioning with Vertical travel of up to 45.9" (116.6 cm); Single-Handed Vertical Adjustment Locking Mechanism
- Articulating swivel arm ($\pm 90^\circ$) allows variety of angled positions
- Receptor Rotation: $\pm 90^\circ$ allows angle flexibility
- Assists in acquiring cross-table lateral exams, as well as many other specialty views
- Receptor Tilt: 0° through 90°
- Quick Locking 5" Caster Wheels

*Specifications are subject to change

www.questxray.com



Improve Your Positioning For
Faster & Easier Examinations